



SkinTyte II™

Deep heating for youthful skin appearance




SCITON®

SkinTyte II™

the speed you need to tyten
large treatment areas

SkinTyte II™ utilizes a range of infrared energy to deeply heat dermal collagen while protecting the entire treated area with sapphire contact cooling. This process promotes contraction of the collagen, initiating the body's natural healing process. Following treatment, a new foundation of collagen is created leading to increased skin firmness.

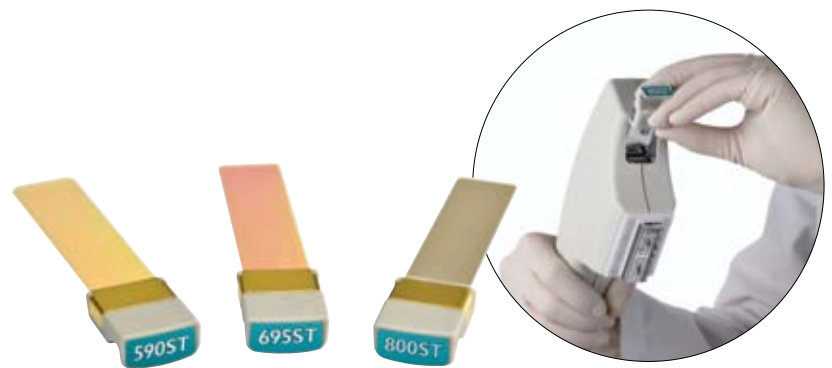

SkinTyte II Smart Filters™



Smoothie Adapters™

SkinTyte II For You:

SkinTyte II delivers the fastest, safest and most economic skin firming treatment solution on the market.

Excellent, Consistent & Predictable Results

Uniform heating targets deeply and effectively to a desired area, immediately stimulating collagen remodeling and contraction.

Fast

With a large 15 x 45 mm spot size, a rapid repetition rate and patient tailored filters, you can treat more patients quickly with SkinTyte II. The magnetically attached Smoothie Adapters facilitate treating with a rapid motion technique.

No Consumables

SkinTyte II does not have any consumables, making it the most cost-effective tissue heating device on the market.

“ A skin-tightening procedure on a large zone like the abdomen can take 30 minutes to an hour with other technologies. We can now get high-quality results in 15-20 minutes with the SkinTyte II. ”

Patrick Bitter, Jr., MD



Before



After 3 Treatments

Photos courtesy of Jody Levine, MD



Before



After 1 Treatment

Photos courtesy of Gloria Lopez, MD

SkinTyte II For Your Patients:

With the SkinTyte II, you can offer your patients a highly requested, effective and comfortable procedure with no downtime.

Powerful Cooling for Maximum Comfort

Precision cooling with SkinTyte II maximizes comfort and safety throughout the treatment. You can cool the skin's surface to precise temperatures from 0 °C to 30 °C.

Safe for All Skin Types

SkinTyte II offers multiple, patient tailored treatment choices for skin types I-VI.

No Downtime

Patients are able to resume their regular activity immediately after treatment making it an extremely popular and in-demand treatment.

The Sciton Edge

Complete Aesthetic Solution – Sciton's SkinTyte II™ comes as part of Sciton's JOULE™ platform or can be included on a standalone BBLs™ BroadBand Light system. JOULE provides the widest array of aesthetic and surgical procedures available on a single, expandable system. The standalone BBLs gives physicians the opportunity to utilize the power of pulsed light to provide a broad range of aesthetic procedures on a single, compact system.

A Secure Investment – JOULE is enabled for tomorrow's technologies and applications so your laser system will not become obsolete. The laser system provides an upgrade path as procedures evolve. Other aesthetic treatment modules can be added to your JOULE, increasing your incremental revenue without the purchase of another full-priced system.

TopFlight™ (U.S. only) – We have partnered with experienced third-party professionals to guide you through the integral aspects of building a cosmetic practice, including business development, marketing strategies, product training, and certification. TopFlight membership is included with most Sciton laser purchases, but please first consult Sciton to confirm your membership.

Practice Support – As part of Sciton's commitment to your success, we provide practice assistance through a Practice Support Kit. This kit includes a comprehensive set of materials that will enable you to easily incorporate the new aesthetic treatments into your practice and to effectively market the treatments to prospective patients.

Clinical In-Service – For non-surgical module purchases or upgrades, Sciton will provide all customers an in-depth clinical training with a credentialed clinical educator. This training is necessary to ensure you and your staff learn proper technique, utilize the most up-to-date Safe-Start protocols, and achieve maximum proficiency when operating your new Sciton system.

Physician Locator (Available in certain countries) – With the purchase of a Sciton laser, you will be listed on Sciton.com's Online Physician Locator. When patients in your area enter their city or zip code into the Physician Locator, they will be given your contact information as being a nearby physician who offers Sciton treatments.

Specifications

Energy Source	Near Infrared – Infrared
Output	Up to 25 W/cm ²
Cooling Method	Continuous thermoelectric cooling – sapphire
Skin Cooling	Adjustable from 0 °C – 30 °C
Spot Size	15 x 45 mm
Smoothie Adapters	15 x 45 mm, 15 x 15 mm



visit us: www.sciton.com

925 Commercial Street, Palo Alto, California 94303
Phone: (888) 646-6999 • Email: info@sciton.com
www.sciton.com



©2011 Sciton, Inc. All rights reserved. Sciton is a registered trademark and SkinTyte II, Smoothie Adapters, and Smart Filters are a trademark of Sciton, Inc.
Printed in USA

2600-029-02 Rev. F