

halo
TRIBRID™

The World's Only
Triple-Wavelength Resurfacing Laser



 SCITON®
BECAUSE RESULTS MATTER



THE EVOLUTION OF SKIN RENEWAL

What If One Treatment Could Do It All?

Imagine delivering dramatic skin renewal—refining texture, correcting pigment, and improving skin quality—all with one handpiece and minimal downtime. Until now, providers have long relied on Sciton's three powerful technologies—HALO®, MOXI™, and Erbium Resurfacing treatments to achieve all three.

Enter HALO® TRIBRID™: a next-generation solution that unites all three *independently* controlled wavelengths—2940 nm, 1927 nm, and 1470 nm—into one powerful treatment on the JOULE®X platform.

Designed to reach key layers of the skin to address the broadest range of skin concerns in a single session, HALO TRIBRID is an evolution in skin renewal—designed to help you deliver better outcomes, in less time, with the confidence only a trusted platform can offer.



The Power of Three: One Handpiece, Unlimited Possibilities

Your patients are unique. And so are their skin concerns.

Until now, providers had to stack multiple treatments to achieve comprehensive results.

HALO TRIBRID changes that.

One device. One session. Three award-winning Sciton technologies—HALO, MOXI, and Erbium Resurfacing.

Together, they give you the power to deliver customized treatments, exceptional results, and the efficiency patients demand, backed by the craftsmanship you trust from Sciton.

2940 nm: Impressive Textural Improvement


This Er:YAG wavelength delivers controlled superficial ablation—**tunable from 20 to 150 microns—for independent, or combination** delivery to refine the appearance of texture, reduce visible imperfections, and accelerate recovery by promoting faster clearance of Microscopic Epidermal Necrotic Debris (MENDs).

1927 nm: Brightening and Pigment Polishing

Now more **powerful with 2x the available adjustable energy**, 10-40* mj, this non-ablative wavelength gently targets the appearance of superficial pigment, enhances skin luminosity, and delivers a refreshed glow that patients notice right away.

1470 nm: Deep Dermal Renewal

Reaching 300-700 microns and capable of **delivering up to 50% faster treatments**, this enhanced, precision-engineered 1470 nm wavelength penetrates the reticular dermis to address the appearance of deeper concerns—**setting the stage for long-term improvement** in tone and texture.

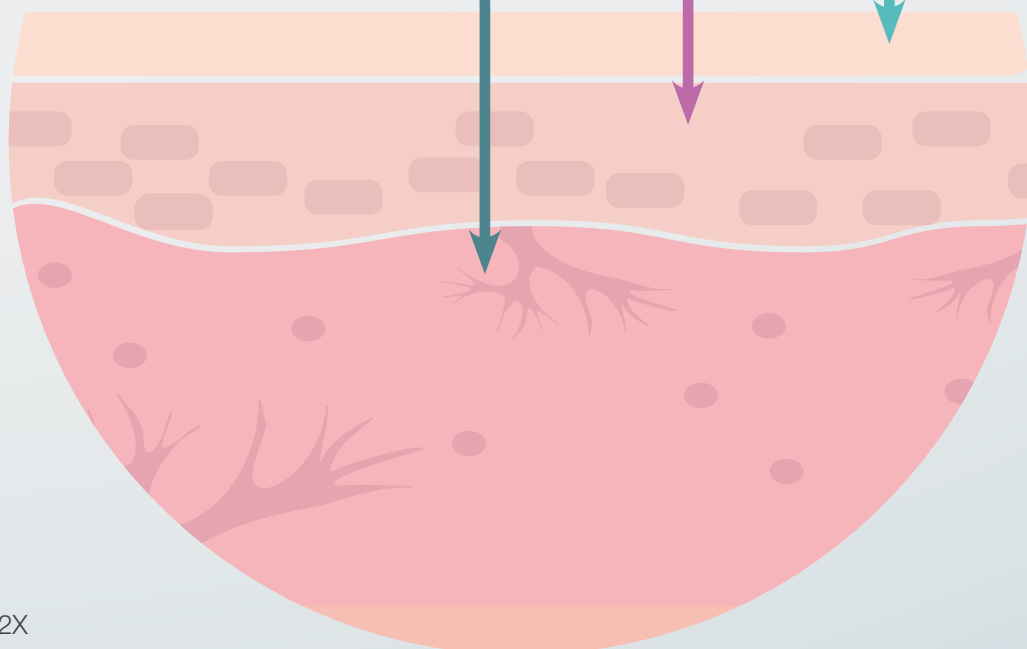

*When 20 mj is delivered 2X



DESIGNED FOR REAL-WORLD DEMANDS

HALO® TRIBRID™ was designed with you in mind: to deliver resurfacing treatments more easily, safely, and efficiently—with results patients will share, remember, and return for. Every detail, from handpiece design to software interface, is designed to help you do exactly that.

One Handpiece. Everything You Need.

The HALO TRIBRID handpiece integrates air cooling and smoke evacuation, so treatments can be performed by a single provider—no assistant required. It's powerful, lightweight, and designed to keep you moving.

Smart Tech That Thinks Ahead

With Dynamic Thermal Optimization (DTO), wavelength automatic adjustments are made in real-time based on the patient's skin temperature—ensuring safe, precise, and consistent energy delivery with every pass.

Easy. Clean. Smart.

The Vortex Adapter, a single-use tip designed specifically for the HALO TRIBRID handpiece, supports proper energy delivery and optimal hygiene, while making setup and cleanup simple.

Real-Time Navigation for Consistent Coverage


The optical navigation system tracks your movement and adjusts scan speed to maintain even coverage—reducing user error and helping you deliver a smooth, uniform result from start to finish.





Software That Works Like You Do

The new user interface is streamlined, intuitive, and designed to keep up with your workflow. Smart navigation, real-time feedback, and effective safe-start settings reduce guesswork and help your team stay efficient—even on busy days.

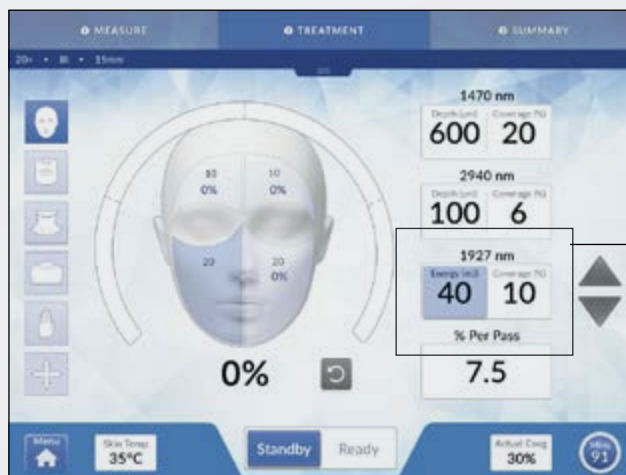

Now more powerful with 2x the available adjustable energy, 10-40* mj.

*When 20 mj is delivered 2x.

2940 nm Erbium is tunable to depths of 100 microns when delivered in combination or 150 microns for independent delivery.



A SMARTER WAY TO WORK. A BETTER WAY TO TREAT

Your time is valuable—and so is your patients'. HALO® TRIBRID™ simplifies your day without compromising results. Less setup, fewer steps, and more impact in every session.

Why Practices Choose HALO® TRIBRID™



One Handpiece, Triple the Impact. Reduce chair time, streamline scheduling, and boost throughput—without extra complexity.



Simplified Workflow. Integrated suction and cooling eliminate the need for additional devices or staff.



Broader Reach, Stronger ROI. Treat all skin types and ages with a single system that delivers year-round value.




Fewer Visits, More Loyalty. Fast, visible results mean quicker rebooking and stronger word-of-mouth.

Why Patient Keep COMING BACK



Three Wavelengths, One Appointment.

Correct, polish, and refine in one efficient session.

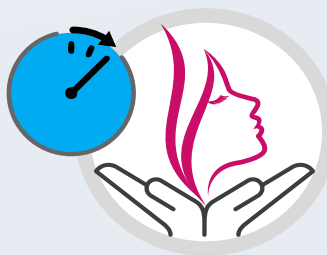

Low Downtime, High Reward.

Most resume normal activity in days, not weeks, with skin that looks and feels younger and healthier.



Safe, Tailored Care.

Adjustable settings for all skin types and goals—so you never compromise on comfort or results.



REAL RESULTS WITH HALO® TRIBRID™

Before & Afters that Build Your Business. Results don't just speak for themselves—they tell your story. With HALO® TRIBRID™, patients notice visible changes from the very first session. That means more trust, faster rebooking, and the kind of before & afters that turn heads—and grow your practice.



HALO® TRIBRID™ | 1 Week Post 1 Treatment



HALO® TRIBRID™ | 1 Month Post 1 Treatment



HALO® TRIBRID™ | 1 Month Post 1 Treatment



HALO® TRIBRID™ | 1 Month Post 1 Treatment



HALO® TRIBRID™ | 1 Week Post 1 Treatment



HALO® TRIBRID™ | 3 Week Post 1 Treatment

Photos Collection Courtesy of
Sherrif F. Ibrahim, MD, PhD | Rochester Dermatologic Surgery



h Post 1 Treatment



HALO® TRIBRID™ | 1 Month Post 1 Treatment



h Post 1 Treatment



HALO® TRIBRID™ | 2 Weeks Post 1 Treatment



s Post 1 Treatment




HALO® TRIBRID™ + BBL HEROic | 1 Month Post 1 Treatment

THE SCITON PARTNERSHIP

At Sciton, we don't just build laser and light-based systems—we build partnerships that last. From the very beginning, our mission has been simple: Do The Right Thing. That commitment has guided us for decades, shaping how we design technology, support providers, and drive innovation in the aesthetics industry.

When you choose Sciton, you're investing in more than just a device. You're choosing best-in-class technology, a future-proof investment, and a team dedicated to your success.

Ready to elevate your practice? Contact us today to learn more about partnering with Sciton.

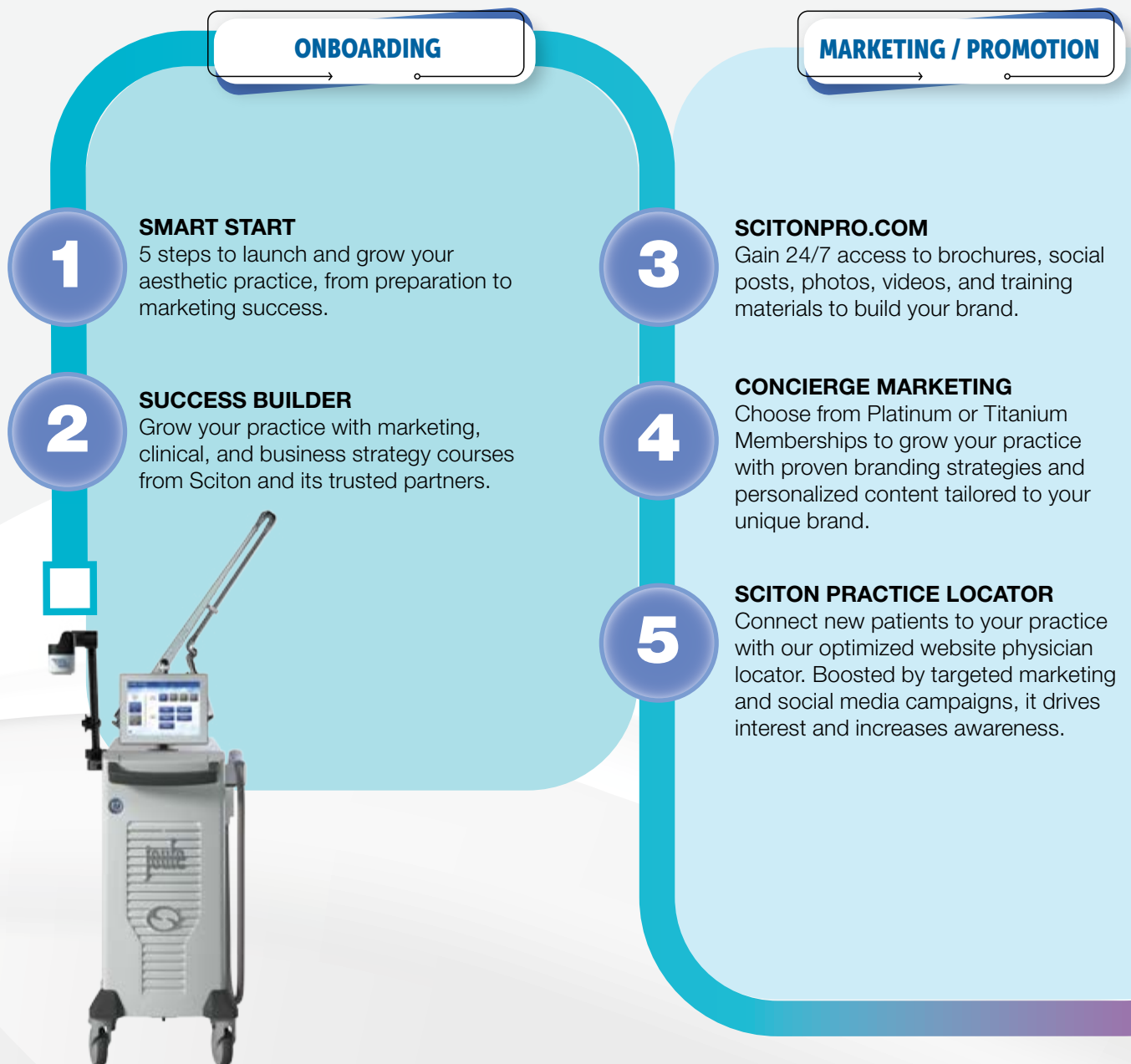

EDUCATION

6

VIDEO LIBRARY

Access quick tip videos right from your platform.

7

CLINICAL EDUCATION

Explore advanced techniques, profitability strategies, and innovative ideas through videos and training materials.

8

EDUCATION AND NETWORKING

Participate in Sciton Supper Clubs, 2-day Skin CEO seminars, and live preceptorships for ongoing professional development.

9

KOL NETWORK

Tap into a network of industry leaders for ongoing support and collaboration to enhance patient outcomes and practice efficiency.

10

WEBINARS

Stay updated on key topics and sharpen your skills with insights from providers and KOLs.

11

SERVICE NETWORK

Keep your platform running at peak performance with support from experienced technicians.

ONGOING SERVICE & SUPPORT

HALO® TRIBRID™ **YOUR COMPETITIVE ADVANTAGE STARTS WITH SCITON**

HALO® TRIBRID™ isn't just another laser treatment—it's the result of decades of innovation, provider feedback, and a relentless pursuit of better outcomes. While other companies outsource manufacturing or compromise on quality, Sciton's award-winning devices on the mJOULE™ and JOULE®X platforms are handcrafted in Palo Alto, CA, with the craftsmanship and quality you can trust.

With HALO TRIBRID, you're not just offering patients a powerful new treatment—you're investing in technology that elevates your outcomes, strengthens your reputation, and drives lasting ROI.

ONE HANDPIECE. THREE WAVELENGTHS. UNLIMITED POSSIBILITIES.

YOU DESERVE
QUALITY
YOU DESERVE
INNOVATION



Corporate Headquarters

925 Commercial Street, Palo Alto, California 94303
Phone: (888) 646-6999 • Email: info@sciton.com
www.sciton.com

Sciton Worldwide

Australia, Canada, Japan, South Korea,
United Kingdom, U.S.A.