mjoule.

The Most Versatile Aesthetic Solution







Affordable quality and versatility in a compact footprint. Meet mJOULE, the award-winning powerhouse platform with treatments for every patient.

BECAUSE

YOU DESERVE

at any budget

SSCITON'

mjoule

MOXI





with in-demand, patient-favorite brands



Clear

Clear Stlk.




SYNERGISTIC TREATMENTS performed alone or together for transformative skin-changing results

mJOULE - The Complete Aesthetic Solution | 3



Whether you're a new or well-established medical practice or med spa, Sciton's mJOULE[™] is the perfect solution to expand your business. Deliver unmatched before and afters for every patient with the platform's consumer-recognized brands and combination treatments.

TARGETED POWERHOUSE MODALITIES

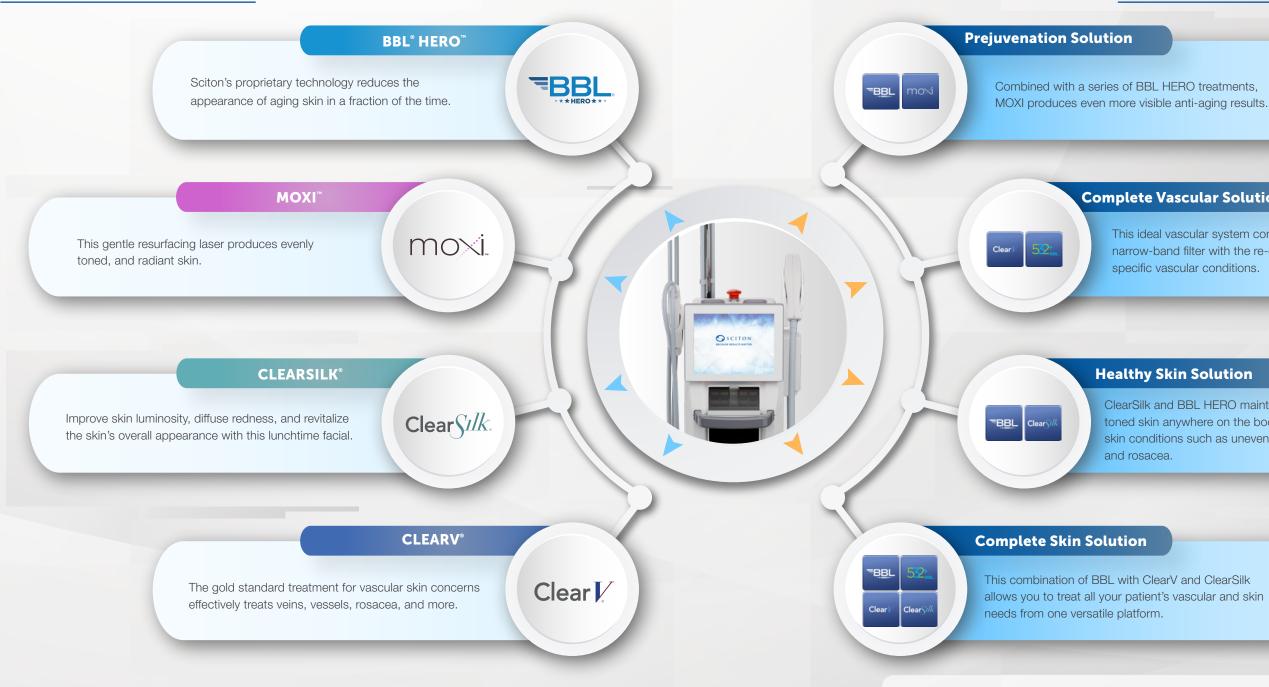

1+1=3 RESULTS – CONFIGURE YOUR POWERHOUSE SYSTEM FOR EXPONENTIAL RESULTS

MOXI produces even more visible anti-aging results.

Complete Vascular Solution

This ideal vascular system combines the new 532+ BBL narrow-band filter with the re-engineered ClearV to treat specific vascular conditions.

Healthy Skin Solution

ClearSilk and BBL HERO maintain healthy and even toned skin anywhere on the body by effectively treating skin conditions such as uneven texture, dull complexion, and rosacea.

Create your dynamic powerhouse system with mJOULE today!



FOREVER CHANGE HOW YOU TREAT SKIN FROM HEAD TO TOE

BBL[®] HERO[™] allows practitioners to successfully treat large areas like the back, legs, and arms in as little as three minutes, thanks to the High Energy Rapid Output (HERO). HERO delivers even better results, allowing for faster treatment and improved patient comfort.

What else makes this device so popular?

With the broadest range of wavelengths on the market, HERO treats pigmented and vascular lesions, acne, and improves the appearance of aged and damaged skin. That's why it continues to be one of the most in-demand treatments today.

Skin Types				
------------	--	--	--	--

INDICATIONS





BBL HERO | Post 3 Treatments Photos Courtesy of Steven Swengel, MD

BBL HERO | Post 1 Treatment

BBL HERO | Post 6 Treatments Photos Courtesy of Sharon Grasso, LE

Dual-Lamp Technology

Dual-lamp technology extends the life of your handpiece and improves performance.



BBL HERO | Post 1 Treatment Photos Courtesy of Jason Pozner, MD

mJOULE - The Complete Aesthetic Solution | 7


A GENTLE APPROACH TO RADIANT SKIN

Patients around the globe love MOXI[™]—the gentle yet powerful 1927nm fractionated laser designed with everyone in mind. It's ideal for those looking for a preventative treatment to help delay the appearance of aging, or those who want to maintain beautiful skin with little to no downtime or discomfort.

In a quick 12-15 minute procedure, MOXI comfortably delivers non-ablative laser energy to correct the initial signs of sun damage and aging—no matter the season or skin type—for a glowing and clear complexion.



INDICATIONS

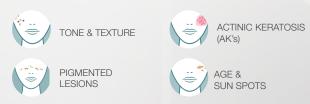

Optical Navigation

The high-speed optical navigation system tracks the motion of the handpiece hundreds of times per second for extreme precision.

Adjustable Beam Placement (ABP)

Users can change the beam scan width, selecting between 15, 7, or 4 mm widths to effectively treat small areas and target lines and scars with incredible precision.

Visual Speed Indicator

The speed indicator lights provide visual feedback to help providers deliver treatments at the optimal speed.



MOXI | Post 2 Treatments Photos Courtesy of Lynn Nieman, MD MOXI | Post 1 Treatment Photos Courtesy of Jason Pozner, MD



The intuitive software allows users to select their treatment preferences from three modes of operation: Manual, Quick Start, and Measure mode.

• Manual: Select your own settings.

MOXI Software

- Quick Start: Choose from three preset treatment levels.
- Measure Mode: Use guided treatment area measurement for even energy delivery.





MOXI | Post 3 Treatments Photos Courtesy of AesthetiCare

Clear Stilk. REVEAL HEALTHY, FLAWLESS SKIN THAT IMPROVES WITH EACH TREATMENT

Your patients are looking for a quick skin treatment with no downtime, but with real results.

ClearSilk[®] is a uniform, non-ablative "lunchtime facial" ideal for first-time laser patients and those who want to maintain a youthful appearance without the recovery time of an ablative resurfacing treatment.

ClearSilk treats redness and rosacea, fine lines and wrinkles, acne scars, and more by depositing heat deep beneath the skin's surface-stimulating fibroblast activity and the body's natural repair process for a healthier, more vibrant complexion. The appearance of skin continues to improve with each treatment.

All Seasons					
Skin Types		IV	V	VI	

INDICATIONS



LESIONS

PIGMENTED

FINE LINES & WRINKLES



Increased Spot Size

Speed up treatment times and uniformly heat the skin with spot sizes between 7 and 10 mm.

Visual Distance Guide

Maintain an appropriate treatment distance from the skin for even heating.

Built-in Temperature Sensing

Get real-time visual feedback from an internal temperature gauge that accurately measures skin temperature.

ClearSilk



ClearSilk | 6 Months Post 3 Treatments Photos Courtesy of Yonick Plastic Surgery

ClearSilk | 3 Treatments Photos Courtesy of Rebecca Gelber, MD

ClearSilk | 6 Months Post 3 Treatments Photos Courtesy of Peter Castillo, MD

Indication Driven Software

The handpiece has special recognition sensors that allow you to directly access the correct application software and safe start treatment settings when taken out of the cradle.

ClearSilk				OSCITON
Fluer Avg. Power (*) 9.6	Steections Skin Types I II III IV V V V	Indications Constraints Constr	3 040 0	Skin Temp 42
«		Standby		

Audible Feedback

Stay alert to how many pulses are delivered in each treatment zone.






ClearSilk | 1 Month Post 4 Treatments Photos Courtesy of Rebecca Gelber, MD

ELEVATE VASCULAR TREATMENTS IN Clear YOUR PRACTICE


At Sciton, we've put a lot of work into developing the ClearV[®] technology to treat vascular skin concerns thoroughly, effectively, and comfortably on any complexion, all year round.

The ClearV treatment eliminates target vessels without damaging the skin. It does so by applying targeted heat to clear away broken capillaries, spider veins, vascular lesions, skin discoloration, and more.

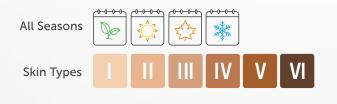

INDICATIONS



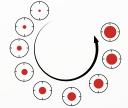

VascuZoom

Instantly change the spot size between 2-10 mm with the adjustable ring directly on the touch display.

Subsurface Vein Illumination (SVI[™])

Enhance visibility and treat hidden subsurface and feeder veins with the push of a button.



Adjustable Integrated Cooling

Ensure patients' comfort and safety with the integrated sapphire cooling tip that optimizes cooling between 5-25 °C before, during, and after treatment. Boost cooling for body areas with a high density of vessels for even more comfort and efficacy.

5-25 °C

Before	During	After

Indication Driven Software

Determine safe start settings and treatment range based on a few key selections. The handpiece has special recognition sensors that allow you to directly access the correct application software when taken out of the cradle.



ClearV | 4 Weeks Post 1 Treatment Photos courtesy of Sherrif Ibrahim, MD, PhD

ClearV | 3 Months Post 1 Treatment Photos courtesy of Sherrif Ibrahim, MD, PhD

ClearV | 8 Weeks Post 2 Treatments Photos Courtesy of Antonio Campo Voegeli, MD





532+ BBL Filter - Only with ClearV

This unique narrow-band filter is the perfect complement to ClearV and in combination with the conical 3 mm adapter, it precisely targets small vascular indications, i.e. telangiectasia, poikiloderma, cherry angiomas and vascular melasma.



3 mm Finesse Adaptor



ClearV | 4 Weeks Post 1 Treatment Photos Courtesy of Antonio Campo Voegeli, MD



ACHIEVE IMPRESSIVE 1+1=3 OUTCOMES WITH **MULTIPLE TREATMENTS**

What if it was possible to achieve significantly superior outcomes for your patients? It is-with the power of ScitonStaX[™] multiple treatments.

StaX[™] with MOXI and BBL HERO effectively treats pigment, tone, and texture in most skin types, including Fitzpatrick type V.



Photos Courtesy of Steven Swengel, MD

Photos Courtesy of Chris Robb, MD, PhD


StaX[™] with ClearSilk and BBL HERO is a powerful solution for reducing rosacea and diffuse redness.



Photos Courtesy of Antonio Campo-Voegeli, MD





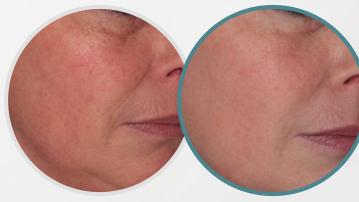

MOXI + BBL HERO | Post 1 Treatment Photos Courtesy of Sherrif Ibrahim, MD, PhD

ClearSilk + BBL HERO | 4 Months Post 4 Treatments Photos Courtesy of Antonio Campo-Voegeli, MD


StaX[™] with ClearV and BBL HERO produces unmatched results for complex vascular conditions such as port-wine stains.




ClearV + BBL HERO | 12 Months Post 4 Treatments Photos Courtesy of Antonio Campo-Voegeli, MD



ClearV + BBL HERO | 6 Months Post 8 Treatments Photos Courtesy of John Cooksey, MD

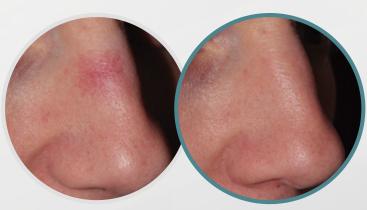


ClearV + BBL HERO | 4 Weeks Post 1 Treatment Photos Courtesy of Antonio Campo Voegeli, MD

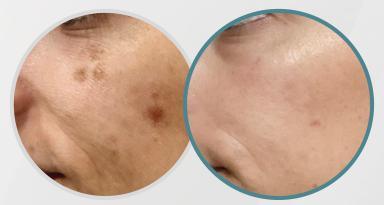

MOXI + BBL HERO | Post 1 Treatments Photos Courtesy of Sharon Grasso, Permanent Touch Cosmetic



ClearSilk + BBL HERO | 1 Month Post 2 Treatments Photos Courtesy of Antonio Campo-Voegeli, MD

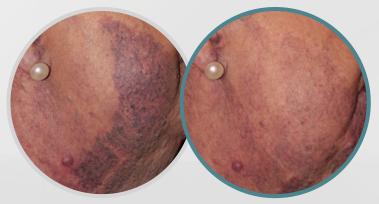

ClearV + BBL HERO | Post 2 Treatment Photos Courtesy of Antonio Campo Voegeli, MD

mJOULE - The Complete Aesthetic Solution | 15

SCITON'S INNOVATIVE TECHNOLOGY IMPROVES LIVES

It's not easy finding a reliable platform that can successfully treat **all your patient's skin concerns**. Especially one that ensures a continuous revenue stream and a solid ROI for your practice. Sciton's mJOULE is a **powerful, revenue-generating workstation**—now more expandable and customizable than ever with world-class modules for all skin types, all year long.

As with all our platforms, the mJOULE is hand-built and designed with exceptional craftsmanship and long-lasting quality—and it's right in line with our mission of developing efficacious and innovative technology that **improves the lives of our customers and patients worldwide.**

YOU DESERVE BETTER | YOU DESERVE SCITON

SPECIFICATIONS

MOXI™

PIQ.	
Wavelength	1927 nm
Delivery	Fiber
Output	1-20 mJ
Power	12 Watt
Treatment Modes	Manual, Quick Start, Measure
Skin Sensor	Optical Tracking
Beam Size	200-300 µm
Treament Depth	Up to 200 µm
Density	Variable

CLEARV® / CLEARSILK®

	Clear V	Clear Silk
Wavelength	1064 nm	1064 nm
Delivery	Fiber	Fiber
Fluence	30-400 J/cm ²	2-27 J/cm ²
Pulse Width	5-200 ms	0.3-1.0 ms
Repetition Rate	1-15 Hz	1-10 Hz
Spot Size	3-10 mm	7-10 mm
Skin Cooling	5-25 °C via Sapphire Tip	Displayed via Handpiece LEDs and on Laser Screen



BBL® HERO™

Wavelength	420 – 1400 nm
Filters	420 nm, 515 nm, 532+ nm*, 560 nm, 590 nm 640 nm, 695
	nm, 800 nm
Output	Up to 30 J/cm ²
Pulse Width	Up to 200 ms
Repetition Rate	1-20 Hz
Cooling Method	Continuous thermoelectric sapphire plates
Skin Cooling	Adjustable from 0 - 30° C in 1 °C increments
Spot Size	15 x 45 mm
Finesse Adapters	15 x 15 mm square, 11 mm, 7 mm & 3 mm round
*Only av	ailable with purchase of ClearV and/or ClearSilk.

Unly available with purchase of cleary and/or clea



Corporate Headquarters



925 Commercial Street, Palo Alto, California 94303 Phone: (888) 646-6999 • Email: info@sciton.com www.sciton.com

Sciton Worldwide

Australia, Canada, Japan, United Kingdom

©2023 Sciton, Inc. All rights reserved. Sciton, BroadBand Light (BBL), ClearSilk, and ClearV are registered trademarks and HERO and mJOULE are trademarks of Sciton, Inc.

