

## The Sciton Edge

- **Complete Aesthetic Solution** – Sciton's JOULE™ laser platform provides the widest array of aesthetic and surgical procedures available on a single, scalable system. ProLipo PLUS™ is one of several laser modules available on JOULE. Other aesthetic treatment modules can be added to your JOULE, increasing your incremental revenue without the purchase of another full-priced system.
- **A Secure Investment** – JOULE is enabled for tomorrow's technologies and applications so your laser system will not become obsolete. The laser system provides an upgrade path as procedures evolve.
- **Practice Support** – As part of Sciton's commitment to your success, we provide practice assistance through a Practice Support Kit. This kit includes a comprehensive set of materials that will enable you to easily incorporate the new aesthetic treatments into your practice and to effectively market the treatments to prospective patients.
- **Surgical Training** – Every surgical module purchase or upgrade includes training with a Sciton-approved physician. Surgical training will cover the use, proper technique, and application of the ProLipo PLUS for laser-assisted lipolysis.
- **Physician Locator** (Available in certain countries) – With the purchase of a Sciton laser, you will be listed on Sciton.com's Online Physician Locator. When patients in your area enter their city or zip code into the Physician Locator, they will be given your contact information as being a nearby physician who offers Sciton treatments.

## Specifications

Wavelength	Nd:YAG 1064 nm and/or 1319 nm
Power (max)	40 W/1064 nm and/or 40 W/1319 nm
Repetition Rate (max)	60 Hz
Pulse Duration	10 - 250 μs

## Accessories

Fiber	600, 1000 micron (Open Source Fibers)
Cannula	14 GA: 200, 250, 300 mm 18 GA: 100, 150, 200 mm
TempASSURE & Cannula	14 GA: 150, 200, 250 mm



TempASSURE™



visit us: [www.sciton.com](http://www.sciton.com)

925 Commercial Street, Palo Alto, California 94303  
Phone: (888) 646-6999 • Email: [info@sciton.com](mailto:info@sciton.com)  
[www.sciton.com](http://www.sciton.com)

©2010 Sciton, Inc. All rights reserved. Sciton is a registered trademark and JOULE, ProLipo PLUS and TempASSURE are trademarks of Sciton, Inc.  
Printed in USA 2600-029-00 Rev. E



## ProLipo PLUS™ + TempASSURE™

High-Power, Laser-Assisted Lipolysis Delivered Safely



# Less FAT. More PATIENTS.

ProLipo *PLUS*™ is a minimally invasive, ultra high-powered laser-assisted lipolysis (LAL) system for melting unwanted fat and coagulating tissue resulting in a tighter, slimmer body. ProLipo *PLUS* employs a dual wavelength capability (1064 nm and 1319 nm), providing the choice to use each wavelength individually or to combine the benefits of both for optimal body sculpting results.

With up to 40 watts of power at 1064 nm and up to 40 watts at 1319 nm, ProLipo *PLUS* is the most powerful and versatile laser-assisted lipolysis system available. Using Sciton's TempASSURE during ProLipo *PLUS* treatment allows this high power energy to be delivered safely and effectively by measuring internal tissue temperature precisely, frequently, and in real-time. This gives physicians a new level of confidence by directly monitoring thermal changes in the adipose tissue in order to identify that the optimal endpoint has been achieved. Further, based on the level of desired tissue heating, it allows the physician to decide how much energy to deliver to a specific treatment area for maximum efficacy.

## ProLipo *PLUS* brings you:

### New Revenue Opportunity

Adding ProLipo *PLUS* to your practice opens your door to a much larger patient population. ProLipo *PLUS* will help you turn more consultations into customers.

### Superior Results

Deliver speed and results with best-in-class power output and dual wavelength capability.

### Unprecedented Flexibility

A wide range of customizable treatment parameters enables body sculpting treatments that are personalized for each patient.

### Maximum ROI

Protect your investment with a fully expandable, modular platform and non-proprietary fibers.

*"The ProLipo PLUS treatment with TempASSURE brings in a different patient base than we typically encounter. A patient will be a lost revenue opportunity if the individual doesn't want surgery in the OR and you don't have an appealing alternative."*

**Marc J. Salzman, MD, FACS**

*"You will be hard-pressed to attain the same degree of skin tightening with more conventional forms of liposuction."*

**Neil C. Goodman, MD**

#### FLANKS




Before

After, Two months post 1 Treatment

Photos courtesy of Marc J. Salzman, MD, FACS

#### WAIST AND THIGHS

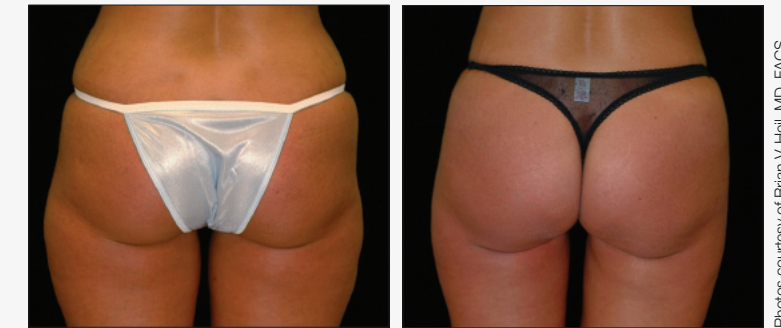


Before

After, Three months post 1 Treatment

Photos courtesy of Brian V. Heil, MD, FACS

#### NECK

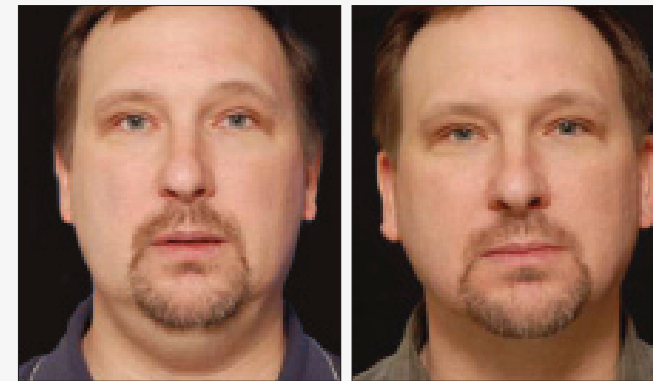


Before

After, Six weeks post 1 Treatment

Photos courtesy of Roger S. Hogue, MD, FVT

#### ARM

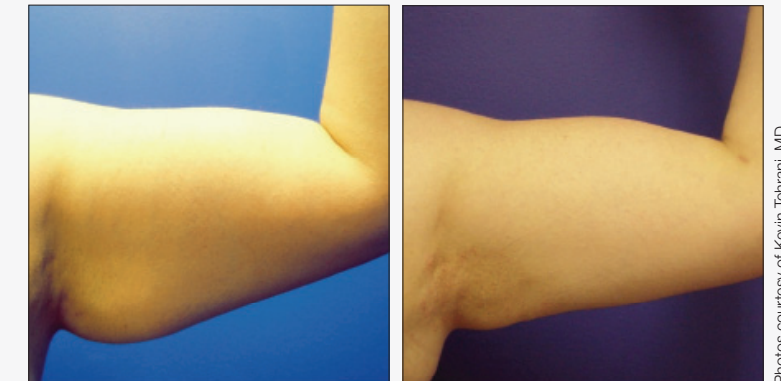

Before

After 1 Treatment

Photos courtesy of Kevin Tehrani, MD

## ProLipo *PLUS* brings your patients:

### Immediately Visible Results

Eliminates fat and tightens skin through heating and coagulation of soft tissue with results that the patient can see immediately.

### Comfortable and Quick Body Contouring

ProLipo treatment with ProLipo *PLUS* is a quick and minimally invasive procedure that allows for a fast recovery time. It does not require general anesthesia, so it can be done in your practice.

### Minimal Downtime

Patients return home the same day and can resume normal activities within a few days.